

# **User Manual**

## **Mandala Breeze**



Please read the manual before using the tools.

Version 2.0

## Table of Contents

<b>1</b>	<b>General Notes .....</b>	<b>2</b>
1.1	Product Feature and Intend Use .....	2
1.2	Warning.....	2
1.3	Warranty .....	2
1.4	Precautions.....	3
1.5	Maintenance.....	3
1.6	Product Includes .....	4
<b>2</b>	<b>Technical Data .....</b>	<b>5</b>
2.1	Specifications.....	5
2.2	Overview and Function .....	6
2.3	Charging Setup.....	6
<b>3.</b>	<b>Fault Diagnosis .....</b>	<b>7</b>



# 1 General Notes

## 1.1 Product Feature and Intend Use

**Thank you for using Light Mandalas Product. Mandala Breeze has been in good condition before leaving factory.**

**Mandala Breeze:** Your personal oasis. Achieve harmony and balance wherever you are. Designed for the rhythm of your life, this sleek, portable device is a masterful blend of technology and wellness. It offers a tranquil escape by utilizing the potent trio of magnetic fields, photon light, and gentle vibration to bring your energy into a state of graceful equilibrium.

Enhance your environment and personal well-being with Mandala Breeze's advanced features. Quantum and Bifilar coils work in unison within this portable device to fine-tune your energy meridians and calibrate your body's natural rhythms. Tailor your experience with our intuitive app and optional Mandala Dynamic Scan, creating a personalized sanctuary of health that adapts to your unique lifestyle needs. With Mandala Breeze, serenity and optimal wellness are at your fingertips, ready to be unfolded wherever life takes you.

## 1.2 Warning

This product is in perfect condition before delivery, please check all accessory when unpacking the unit.

Pay close attention to each safety REMARK and WARNING statement in the manual. This manual contains important operation information.

For your eyes safety, please wear light protective goggles and avoid light to expose to your eyes.

## 1.3 Warranty

Each single product passes through a quality process. Nevertheless, if failures occur, contact the customer service at once. The product is under 1 year warranty effectively according to the delivery date stated on the warranty card.

Warranty does not apply to damages, which result from misuse or neglect. The warranty also expires if you dismantle the product or failed to provide the warranty card. The manufacturer is not liable for consequential damage or in case of a non-intended use of the product.

If a failure occurs during the warranty period, the product will be replaced or repaired. The manufacturer reserves the right to exchange components of the product instead of repairing it. If the failure results from misuse or neglect, the user must pay for all cost that may occur. The freight cost and other costs will be depending on which Term & Conditions criteria the product falls into. **Import tax and Customs clearance duty cost is not included in the term of warranty.** In that case, you may ask for an estimated cost beforehand.

## 1.4 Precautions

1. Please check the device for damage during shipping. Do not use the device if any damage is discovered. If the device is damaged, please contact the dealer or Light Mandalas immediately.
2. Please Confirm the product you have received are in correct voltage and frequency for your country.
3. Mandala Violet Shield needs to be charged at least 2 hours prior to first use.
4. Do not disassemble or modify the device in any case.
5. Do not use the device outdoors when it rains.
6. Place the device in a well-ventilated area.
7. Keep it dry and avoid excessive humidity, heat, or dust.
8. Do not place the device closer than 2 feet from combustible objects.

## 1.5 Maintenance

1. If the device has not been used for a long time, please charge for 2 hours before use.
2. Clean with soft cloth, humid tissue moistened with water or a water-based glass cleaner.
3. To guarantee the fixture works properly, it is essential to keep it clean.

## 1.6 Product Includes




1. Mandala Breeze
2. USB Type-C cable
3. Mandala Breeze's box

## 2 Technical Data

### 2.1 Specifications

Product name:	Mandala Breeze
Color:	White
Material:	ABS plastic + Acrylic plastic
Input Voltage:	DC 5V / 2A
Charging Interface:	Type C USB
Battery Type:	3.7v 500 mAh LiPo Battery
LED Output:	LED RGB 3528
Coil Output:	Bifilar / Quantum
frequency:	0 – 2MHz
Bluetooth:	BLE JDY-23
Length:	45 mm
Width:	45 mm
Height:	18 mm
Weight:	32 g

## 2.2 Overview and Function

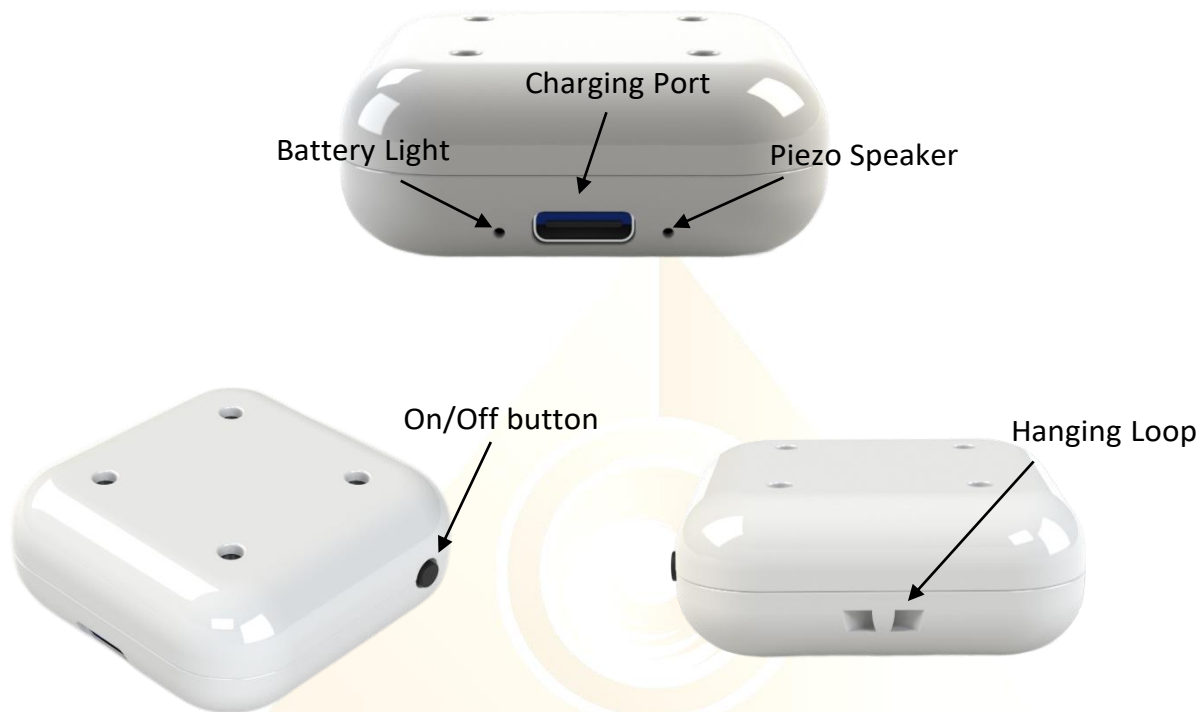

**\*ON/OFF button:** Press once to turn on device and double tap to turn off device.

## 2.3 Charging Setup

1. To charge Mandala Breeze, plug in the **USB Type-C Cable** to the device charging port and connect to **Charging Adapter** and to the power supply. The battery light and device will be lighting up while charging the battery. The battery light will stop blinking when the battery is completely charged. Then we can turn off the device first before turn on device to apply.



- The charging process takes about 2 hours if starting with a fully drained battery and provides over 6 hours of battery life.



*To protect the battery life span, please don't overcharge the battery.*

### 3. Fault Diagnosis

Failures	Remedies
Device cannot turn on	<ul style="list-style-type: none"><li>The battery life is low. Please charge the device.</li><li>Plug in the charger. If it can be used during charging, that means the device needs to be charged.</li></ul>



## Disclaimer

All content found on the Light-mandalas.com Website, Facebook Page, YouTube, including: text, images, video, audio, or other formats are intended for informational, educational and research purposes only and neither Light Mandalas Co. Ltd. nor the publishers of this content take responsibility for possible health consequences of any person or persons reading or following the information in this educational content and is not intended to be a substitute for professional medical advice, diagnosis, treatment or cure. Always seek the advice of your physician or other qualified health professionals with any questions you may have regarding a medical condition.

**“No claims other than stress detection and stress reduction may be made”**

Light Mandalas products are to be used as a Heart Rate Variability Biofeedback system and Light, Color and Sound therapy, designed for stress detection and stress reduction only. The devices do not evaluate any issue other than stress. Stress can come from many sources; this system uses light and color and sound therapies to reduce stress. Remember, there has never been any reported side effects using light, color, and sound therapy anywhere in the world and throughout time.

### Notes on the practices of low-level laser therapy (LLLT) and Biofeedback

Low-level laser therapy (LLLT) is considered investigational by the FDA, and no medical claims of cures are permitted in the US; at this time however, low-level lasers are defined within the scope of complimentary medical practices such as physiotherapy, aesthetic, Traditional Chinese Medicine, and acupuncture, etc.

The FDA does not regulate Bio resonance tools, but they do regulate biofeedback tools and some Bio resonance tools are registered under this category.

\*It is a legal requirement at this time that the user does not make any medical claims of diagnosis or cure.

### Precaution

Light Mandalas products are working on energetic level. The standard benefit is energy clearance. During the process, your body and mind can be uncomfortable. The symptoms associated with overly rapid clearance of toxins. If those detox symptoms are too strong, you may need to stop using the products. You have to allow the clearance or detox to proceed, and you can start to try the products afterwards. However, if the issues still remain the same, it may need to be completely discontinued for your own safety

\*Although UVC is considered the safest of the UV wavelengths please do not place stationary on the bare skin for extended period of time (more than 1 minute) and avoid looking directly into the light source.

**Shall there be any inquiry, please contact us for further support.**

**Light Mandalas Co., Ltd. (Thailand Headquarters):**

- **Light Mandalas Website:** [www.light-mandalas.com](http://www.light-mandalas.com)
- **Official:** [info@light-mandalas.com](mailto:info@light-mandalas.com)
- **Instagram:** [@mandalawave9](https://www.instagram.com/mandalawave9)
- **YouTube Channel:** Light Mandalas Official
- **Technical Support Email:** [cservice@light-mandalas.com](mailto:cservice@light-mandalas.com)

**Local agents contacts:**

- **Japan:** [japan@light-mandalas.com](mailto:japan@light-mandalas.com)
- **Taiwan:** [lovelight777888@gmail.com](mailto:lovelight777888@gmail.com)
- **Indonesia:** [lightseed2020bali@gmail.com](mailto:lightseed2020bali@gmail.com)
- **West Malaysia:** [lightworker@healingtech.org](mailto:lightworker@healingtech.org)
- **East Malaysia:** [malaysia@light-mandalas.com](mailto:malaysia@light-mandalas.com)
- **Hungary:** [horesnyibela@light-mandalas.com](mailto:horesnyibela@light-mandalas.com)
- **Switzerland, Germany, Austria, and the Netherlands:** [info@light-mandalas.ch](mailto:info@light-mandalas.ch)
- **France, Spain:** [france-spain@light-mandalas.com](mailto:france-spain@light-mandalas.com)
- **USA:** [usa@light-mandalas.com](mailto:usa@light-mandalas.com)
- **New Zealand, Australia:** [tachyonhealing.nz@gmail.com](mailto:tachyonhealing.nz@gmail.com)